

Human CXCL10/IP-10 ELISA-KIT

Ref.: hIP-10-EIA-5-KIT

MabTag's ELISA for human Chemokine (C-X-C Motif) Ligand 10 (CXCL10) or Interferon gamma-induced protein 10 (IP-10) contains appropriate reagents sufficient for processing of 5 microplates (5 x 96 wells; 100 µl/well)

For research only. Not for use in diagnostic or therapeutic procedures.

Specificity: human Chemokine (C-X-C Motif) Ligand 10 (CXCL10)

Typical standard curve range: 8 – 500 pg/ml

Detection limit: 2.6 pg/ml

Samples: Culture supernatants, serum, plasma and other body fluids.


For serum and plasma a dilution of $\geq 1:10$ is recommended.



MabTag
GmbH

Content	Working dilution	Storage
1 x vial 500 µl liquid anti-hCXCL10 Capture-Antibody (red cap)	1:100	-20°C
1 x vial 500 µl liquid anti-hCXCL10 Detection-Antibody (yellow cap)	1:100	-20°C
1 x vial 50 ng lyophilized rhCXCL10 Standard (white cap)	customer specific	-20°C
1 x vial 50 µl Poly-HRP-Streptavidin (blue or green cap)	1:1000	-20°C
5 x 96well-Microplate	-	room temperature
1 x 60 ml Coating-Buffer	ready-to-use	2-8°C
3 x 100 ml Blocking-Buffer / Reagent-Diluent	ready-to-use	2-8°C
3 x 100 ml Wash-Buffer (10x concentrated)	100 ml must be filled up to 1L with aqua dest.	2-8°C
1 x 4 ml TMB-Solution A	10 ml TMB-B + 0.5 ml TMB-A	2-8°C
1 x 60 ml TMB-Solution B		protect from light
1 x 30 ml Stop-solution ACID → wear gloves & protective glass!	ready-to-use	2-8°C

Typical standard curve

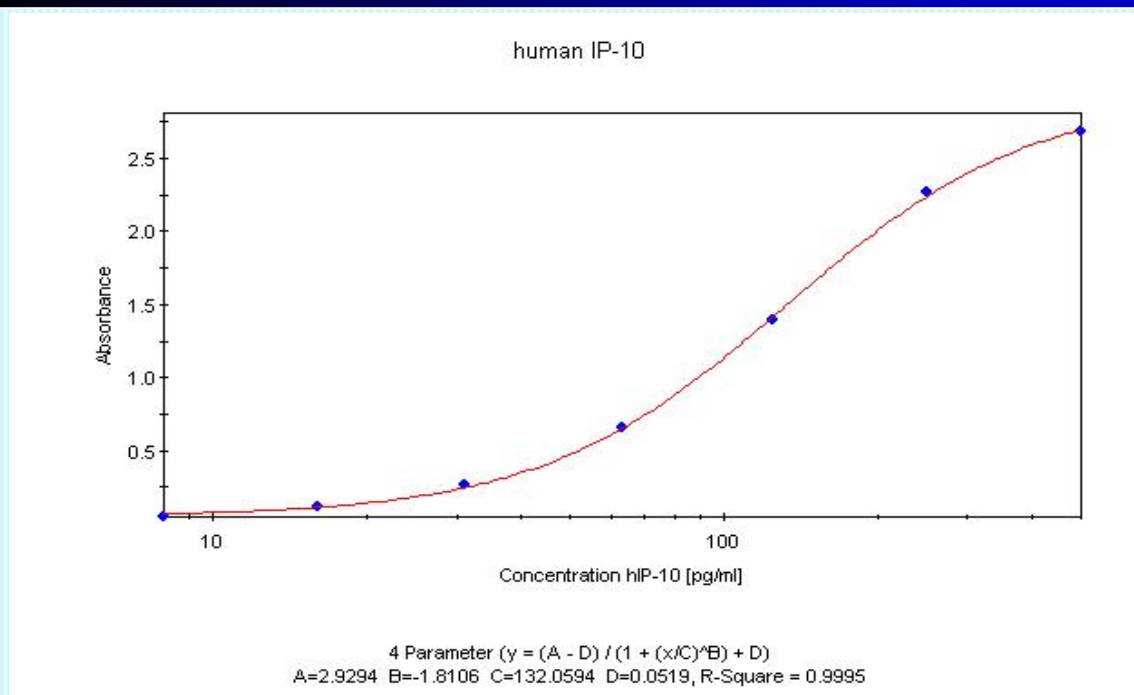

!Spin down all vials before use!

Step	Incubation	Procedure
Coating Capture-antibody	≥ OVERNIGHT at room temperature	Dilute capture-antibody 1:100 in COATING-BUFFER (100 µl capture-antibody in 10 ml COATING-BUFFER). Subsequently transfer 100 µl of this working-solution to each well and incubate.
Remove capture-antibody completely by inverting the microplate and blotting it vigorously against clean paper towels.		
Blocking	≥ 1 Hour at room temperature	Add 300 µl BLOCKING-BUFFER to each well and incubate.
Remove BLOCKING-BUFFER completely by inverting the microplate and blotting it vigorously against clean paper towels.		
Standard & Sample	≥ 2 Hours at room temperature	Dilute standard & samples in REAGENT-DILUENT and transfer 100 µl in the respective wells in duplicates. Standard: Make serial dilutions in REAGENT-DILUENT and begin with a high standard and end with blanks. The standard vial of this set contains 50 ng lyophilized standard . Reconstitute this in exactly 1 ml REAGENT-DILUENT (stock solution = 50 ng/ml) and choose a sufficient high standard working solution for your assay (e.g. prepare a 1:25 dilution for a standard curve beginning with 2000 pg/ml).
Wash 5x vigorously with WASHING-BUFFER and remove resting buffer completely by inverting the microplate and blotting it vigorously against clean paper towels.		
Detection- antibody	≥ 2 Hours at room temperature	Dilute detection-antibody 1:100 in REAGENT-DILUENT (100 µl detection-antibody in 10 ml REAGENT-DILUENT). Subsequently transfer 100 µl of this working-solution to each well and incubate.
Wash 5x vigorously with WASHING-BUFFER and remove resting buffer completely by inverting the microplate and blotting it vigorously against clean paper towels.		
Poly-HRP- Streptavidin	20-30 Min at room temperature	Dilute Poly-HRP-Streptavidin 1:1000 in REAGENT-DILUENT (10 µl in 10 ml REAGENT-DILUENT). Subsequently transfer 100 µl of this working-solution to each well and incubate.
Wash 5x vigorously with WASHING-BUFFER and remove resting buffer completely by inverting the microplate and blotting it vigorously against clean paper towels.		
Substrate solution	Up to 60 Min* at room temperature in the dark	Optionally warm the solution to room temperature before use. Add 100 µl of the SUBSTRATE-SOLUTION to each well and incubate. Control the development of the colour reaction continuously and stop at an appropriate time point.
Stop solution	-	When the enzymatic colour reaction is sufficiently proceeded stop the reaction by adding of 50 µl stop solution. Read the microplate at the substrate-depending wavelength. (e.g. 450 nm for TMB substrate) (if wavelength correction is available, set to 540 nm, 570 nm or 630 nm as reference)

*The speed of enzymatic colour development is influenced by many customer-specific factors. Therefore the incubation time is variable und specific for each test system.

Note:

All incubation steps except TMB substrate could be optionally carried out over-night. Do not use sodium azide-containing solutions, nor add sodium azide to the supplied reagents. Sodium azide inactivates the peroxidase.

Storage:

Specific storage conditions in the table above.

Reconstituted reagents should be stored at -20°C. Please prevent repeated freeze- thaw cycles. Stable for up to 6 months after opening when stored at -20° C. The performance of the unopened reagents is guaranteed until one year after point of delivery.

Precautions for use:

!The stop solution is an acid solution. TMB-Solution A contain H₂O₂ and tetramethylbenzidine (TMB). All Buffers and liquid antibody solutions contain 0.045% (v/v) Proclin®950 as preservative. All these compounds are harmful and cause respiratory, skin and eye irritation. Do not swallow any components of the set/kit (R22). Wear face, eye and hand clothing protection when using this material (S36). Keep out of reach of children (S2). Keep away from food, drink and animal feeding stuff (S13). !These reagents are offered for research purposes only! For *in vitro* use only. MabTag will not be held responsible for patent infringement or other violations that may occur with the use of our products.

ProClin is a registered trade mark of Rohm and Haas Company.